
waterpik™

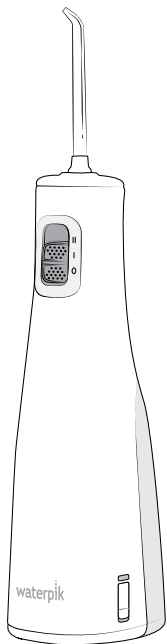Waterpik™ Water Flosser Model WF-20 Series

Idropulsore Waterpik™

Modello Serie WF-20

Irrigador bucal Waterpik™

Modelo Serie WF-20

Waterpik™ Water Flosser

Model WF-20 serie

www.waterpik.com/intl

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place or drop the product into water or other liquid.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Do not handle the user supplied USB power supply (power adapter) with wet hands.
- Unplug the user supplied USB power supply (power adapter) when not in use.

- Check the charging cord and user supplied USB power supply (power adapter) for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only for its intended use as described in these instructions or as recommended by your dental professional.

- Never operate this product if it has a damaged cord or external power supply, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact your local distributor which can be found at www.waterpik.com/intl.
- A user supplied USB-A wall charger (power adapter), safety certified (IPX4, Class II) with a rated output of 5.0 V $\overline{=}$, 1.0 A Max is required to charge this product.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Keep the charging cord away from heated surfaces.
- Never drop or insert any object into any hose or opening.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Fill reservoir with water or other dental professional recommended solutions only.
- Only use tips and accessories recommended by Water Pik, Inc. as identified in this instruction manual.

- If the charging cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product.
- Do not direct water from the water flosser under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.

- Do not use the water flosser for more than five minutes in each two-hour period.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
- There are no user serviceable electrical parts in this product and it does not require maintenance.
- Remove the charging cord from the unit and replace charging port cover prior to filling the reservoir.
- Ensure contacts on the cord, the unit and the power adapter are clean and dry while the unit is being charged.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet. Contact your doctor with any questions.
- Do not expose the battery or appliance to excessive temperatures.

SAVE THESE INSTRUCTIONS

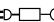
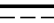
INTENDED USE

Waterpik™ water flossers with PRECISIONPULSE TECHNOLOGY are proven to deliver an effective combination of water pressure and pulsations to remove plaque bacteria that can cause gingivitis, cavities, and bad breath.


Product Description	7
Getting Started	8
Cleaning and Troubleshooting Your Waterpik™ Water Flosser	11
Battery Care, Battery Disposal/Recycling and Battery Removal	12
Warranty Information	13

Symbols Glossary (Not all symbols apply to all products)


Explanation of symbols used on Water Pik, Inc. product labels and packaging. **Applicable symbols for a specific product are found on that product's labels and packaging.**

Symbol	Definition	Symbol	Definition	Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Detachable Power Supply		Do not use this device in a shower, bathtub and do not immerse in any liquid.		Separate collection for waste electrical and electronic equipment.
	Manufacturer		Direct Current (DC) Voltage		Product contains a lithium ion battery		Contains a recyclable Li-ion battery
IPX7 or IPX4	Ingress Protection Code		Alternating Current (AC) Voltage		CE marking of conformity		Contains a recyclable Ni-MH battery
	IEC Class II Equipment		Do not immerse this device in any liquid.		Regulatory Compliance Mark		

1. Power handle
2. Power/pressure switch
3. Battery charge indicator
4. Charge port/cover
5. Tip
6. Tip eject button
7. Reservoir fill door
8. Reservoir

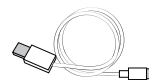
Charging Cable
20033894-2 White

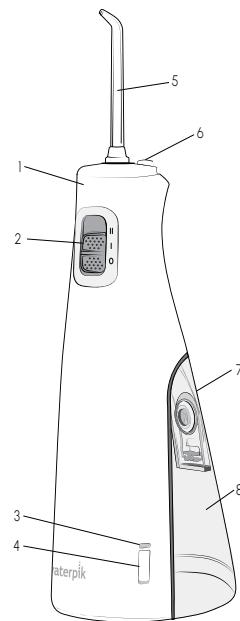
TIPS

- | | |
|--|---|
| A. Orthodontic Tip* (OD-100E)
- Braces/General use | D. Handheld Cordless Precision Tip* (JT-450E)
- General use |
| B. Pik Pocket™ Tip* (PP-100E)
- Periodontal pockets/Furcations | E. Tongue Cleaner** (TC-100E)
- Fresher breath |
| C. Plaque Seeker™ Tip* (PS-100E)
- Implants/Crowns/Bridges/Retainers/General use | F. Travel Plug
- Keeps excess water from draining out while traveling |

NOT ALL TIPS INCLUDED IN ALL MODELS

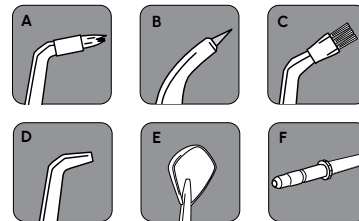*Replace Every 3 months **Replace Every 6 months
Replacement tips may be purchased from your local distributor which can be found at www.waterpik.com/intl.

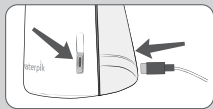Charging Unit Prior to First Use

Connect the charging cable to a USB power adapter (not included), then plug into a wall outlet. Attach the other end of the cable to the unit.

Charge for 4 hours before first use. A full charge lasts up to 9 weeks with typical use.

Replace the charge port cover after charging and ensure it's flush with the handle.

Note: Refer to the IMPORTANT SAFEGUARDS section for required USB power adapter specifications.

Charge Lights

- The LED will flash white to indicate that the unit is approaching a low battery level.
- The LED will turn red when the battery is depleted, indicating that the unit needs to be charged.
- The LED will glow when charging.
- The LED stays lit when fully charged until the charging cable is unplugged from the device.

Filling the Reservoir





Before filling, disconnect charging cable and replace charging port cover.

Lift flip top on reservoir and fill with lukewarm water.

Inserting and Removing Tips

Insert tip into the center of the knob at the top of the water flosser handle. Press firmly until the tip clicks into place. To remove the tip from the handle, press the tip eject button and pull tip from the handle.

CAUTION: DO NOT PRESS TIP EJECT BUTTON WHILE UNIT IS RUNNING.

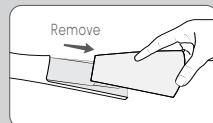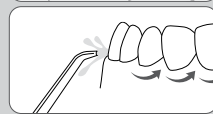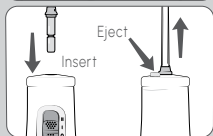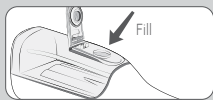Adjusting the Pressure Setting

You can adjust the pressure setting on your water flosser by simply sliding the pressure/power switch on the handle. Slide power switch up 1 click for Low (I) pressure and 2 clicks for High (II) pressure.

Start with the lowest setting and increase to your liking.

Recommended Technique


Lean low over sink and place tip in mouth. Aim the tip toward gumline, and turn the unit ON (I). For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth. Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

Tip:

- LOOK INTO SINK, not in mirror, to reduce splashing.

When Finished

Turn the unit OFF (O). Empty any liquid left in the reservoir. If desired, the reservoir may be removed from the power handle by simply pulling down on the reservoir to remove.

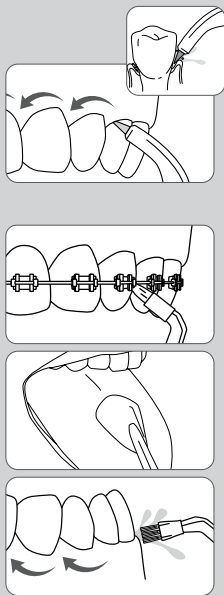NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

SPECIALTY TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets. To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on and continue tracing along the gumline.

Orthodontic Tip


To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

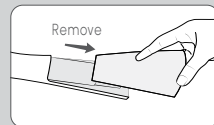Using Mouthwash and Other Solutions

Your Waterpik™ water flosser can be used to deliver mouthwash and antibacterial solutions. For ease of adding other solutions remove the reservoir (see "When Finished" section) and add from fill door on reservoir. If you are adding additional water you may replace the reservoir and fill reservoir fully with the fill door. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

The water reservoir is removable for easy cleaning and is top rack dishwasher safe. We recommend air drying instead of a heated dry cycle.

It is recommended that the reservoir is removed and cleaned weekly.

Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

This process should be done every 1 to 3 months to ensure optimal performance.

Battery Care

If the water flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage. Do not expose the battery or appliance to excessive temperatures.

Battery Disposal/Recycling

This product contains a **rechargeable lithium ion battery**. The battery is not replaceable. At the end of product life, do not dispose of in household garbage or recycling bins. Please recycle at your nearest recycling center or take to hazardous waste collection point to comply with local and state regulations.

Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire. Please contact your local waste administrative department for any questions.

Note: The rechargeable battery is not replaceable, remove only at end of product life. This process is not reversible and voids the product warranty.

Battery Removal



To remove the rechargeable battery for recycling, please discharge the battery fully before removing. Observe basic safety precautions when you follow the procedure below. Be sure to protect your eyes, hands, fingers and surface on which you work.

1. Remove reservoir and wrap unit in towel. Hit with hammer until lower cover comes off.
2. Pull out battery, cut wires. Cut battery wires one at a time.

WARNING:

DO NOT cut both battery wires at the same time.

DO NOT disassemble, puncture, hit, squeeze or throw the battery into fire.

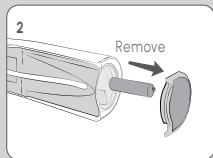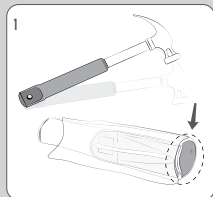Service Maintenance

Waterpik™ water flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all parts and accessories, contact your local distributor at www.waterpik.com/intl. Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit. Still have questions? Contact your local distributor which can be found at www.waterpik.com/intl.

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase (Spain Only: 3 years). **Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim.** Water Pik, Inc. will replace the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories

or consumable parts such as tips, etc. **Use with improper voltage will damage product and VOID warranty.** This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

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©2025 Water Pik, Inc.
**Water Pik, Inc., a subsidiary of
Church & Dwight Co., Inc.**
SOFIBEL SAS
Church & Dwight UK Limited,