


QUICK START GUIDE

Model WF-06


Register your product and receive over £45 in savings.
Waterpik.co.uk/registration

Visit Waterpik.co.uk/WF06
to watch How To Video


1 Unit Prep (First time use only)

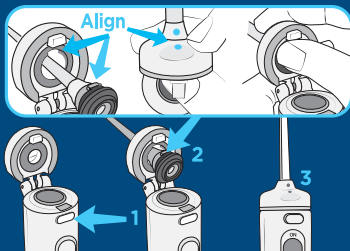

A. Press button to open flip top.¹ Insert tip through top, align notches on tip and handle (dots on tip and flip top will align), press firmly to click tip into place.² Close flip top.³

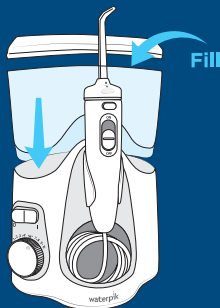

B. Fill reservoir with warm water and replace it firmly on the base.

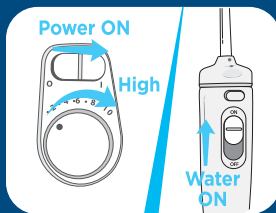


To Prime Unit

C. Turn pressure to high, point tip into sink, turn power and water switches ON until water flows. Turn power and water OFF.

Register your product and receive over £45 in savings.
Waterpik.co.uk/registration

Visit Waterpik.co.uk/WF06
to watch How To Video

2 Water Floss

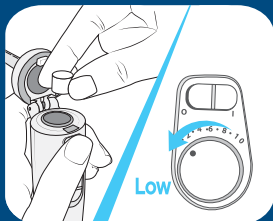


A. Load tablet into handle. Start on low pressure and increase to your liking.



B. With water switch in OFF position, place tip in mouth. Turn power ON.



C. Lean over sink. Turn water switch to ON position.

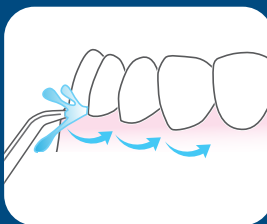

D. Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.