


# waterpik® waterflosser®

## QUICK START GUIDE


Review full instruction book prior to first use

### 1 FIRST USE

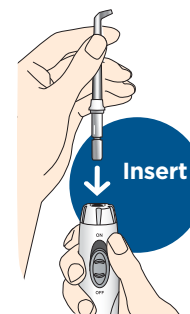

Insert tip firmly until it clicks into place.

Turn pressure to high, point tip into sink, turn power and water switches on until water flows. Turn off.

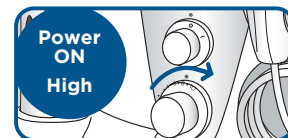

Model WP-950

This device is compatible with global voltages of 100-240VAC, 50/60Hz.

©2016 Water Pik, Inc.

# waterpik® waterflosser®


## QUICK START GUIDE

Review full instruction book prior to first use

### 2 START WATER FLOSSING


Turn power and water on. Start on low pressure and increase to your liking.

**TO AVOID MESS:  
PLACE THE TIP  
IN YOUR MOUTH  
BEFORE TURNING  
UNIT ON**

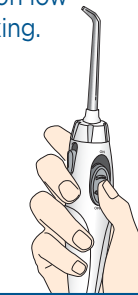


### 3 LEAN OVER SINK

Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.

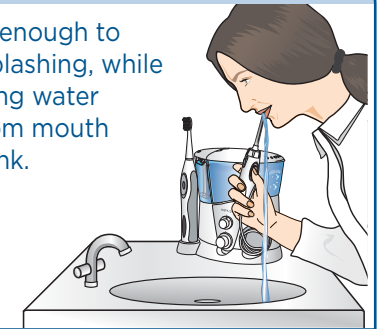

### 4 TECHNIQUE

Aim water just above the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

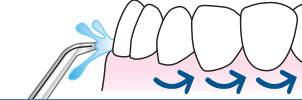

### 5 SEE HOW TO VIDEO: [waterpik.com/HowTo950](http://waterpik.com/HowTo950)

Receive up to £45 in savings on the purchase of additional Waterpik® products.  
[www.waterpik.co.uk/registration](http://www.waterpik.co.uk/registration)

20021915-F AA  
+44 (0) 333 123 5677