


waterpik®

# Waterpik® Water Flosser Model WP-120

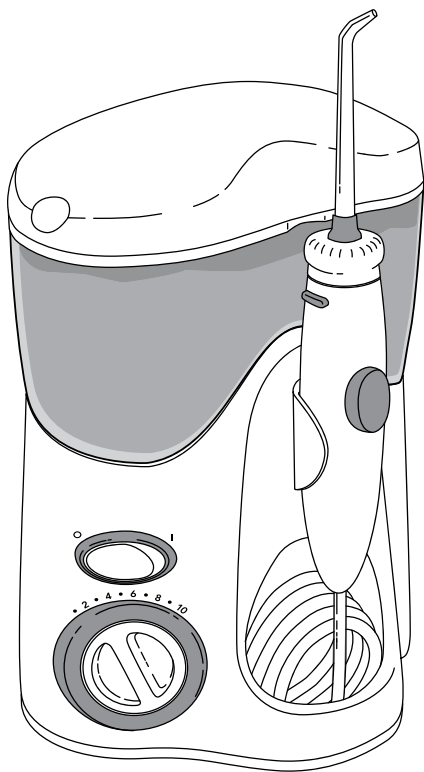

.....  
[www.waterpik.co.uk](http://www.waterpik.co.uk)

## TABLE OF CONTENTS

|   |   |
|---|---|
| Important Safeguards                                      | 3 |
| Product Description                                       | 5 |
| Getting Started   | 6 |
| Cleaning and Troubleshooting Your Waterpik® Water Flosser | 8 |
| Limited Three-Year Warranty                               | 9 |

## IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### READ ALL INSTRUCTIONS BEFORE USING.

### DANGER:

#### To reduce the risk of electrocution:

- Always unplug product after using.
- Do not handle plug with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a bath or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.


### WARNING:

#### To reduce the risk of burns, electrocution, fire, or injury to persons:





- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Waterpik International, Inc. at +44 (0) 844 335 6354.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

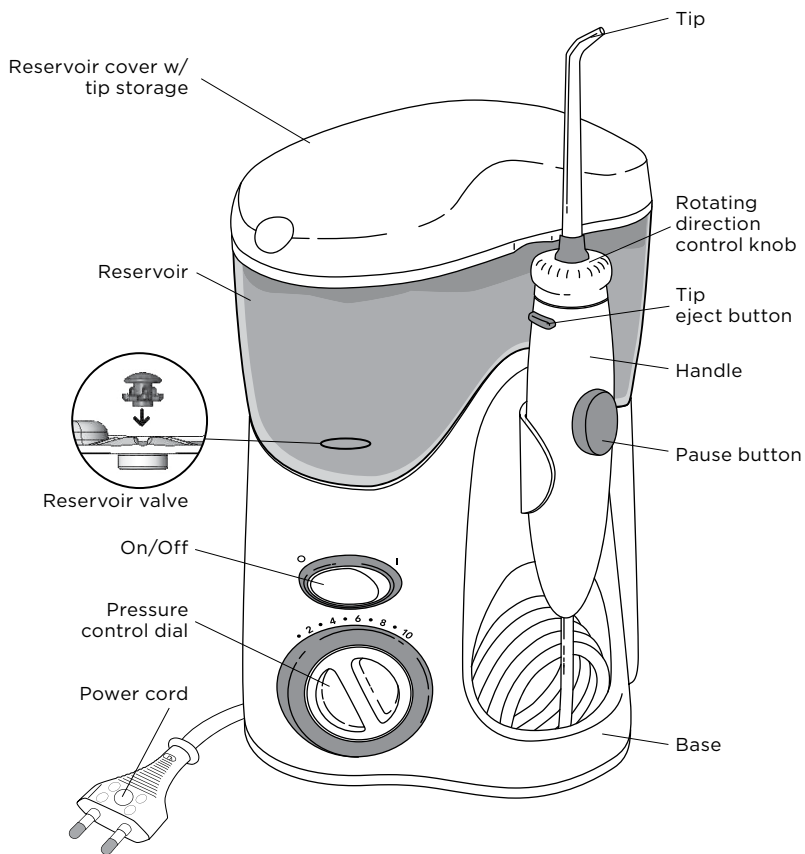- Do not use iodine, bleach, or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
- Do not use this product for more than 5 minutes in each two-hour period.
- Fill reservoir with warm water or other dental professional recommended solution only.
- Children should be supervised to ensure that they do not play with the appliance.
- Closely instruct and supervise children 6 years and above and individuals with special needs in the proper use of this product.

**SAVE THESE INSTRUCTIONS.**



## GETTING STARTED

| Replace Every 3 Months  |  |
|---|--|
| <b>Orthodontic Tip</b>  | <b>BEST FOR</b>  |
|  | <input checked="" type="checkbox"/> Braces<br><input checked="" type="checkbox"/> General Use  |
| <b>Pik Pocket™ Tip</b>  | <b>BEST FOR</b>  |
|  | <input checked="" type="checkbox"/> Periodontal Pockets<br><input checked="" type="checkbox"/> Furcations  |
| <b>Plaque Seeker® Tip</b>   | <b>BEST FOR</b>  |
|  | <input checked="" type="checkbox"/> Implants <input checked="" type="checkbox"/> Bridges<br><input checked="" type="checkbox"/> Crowns <input checked="" type="checkbox"/> General Use |
| <b>Toothbrush Tip</b>   | <b>BEST FOR</b>  |
|  | <input checked="" type="checkbox"/> General Use  |

| Replace Every 6 Months  |  |
|---|--|
| <b>Classic Jet Tip</b>  | <b>BEST FOR</b>                                    |
|  | <input checked="" type="checkbox"/> General Use    |
| <b>Tongue Cleaner</b>   | <b>BEST FOR</b>                                    |
|  | <input checked="" type="checkbox"/> Fresher Breath |

**NOTE:** All tips not included in all models.

Replacement tips/attachments may be purchased directly from Waterpik International, Inc. +44 (0) 844 335 6354 or from your local Water Pik stockist.

### Powering the Unit

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.

### Preparation and Use

Refer to the Quick Start Guide for proper use instructions.

### When Finished

Turn the unit OFF (O). Empty any liquid left in the reservoir. Unplug the unit.


**NOTE:** If you used mouthwash or other solution besides water, see “Using Mouthwash or Other Solutions” section.

## TIP USAGE

### Pik Pocket™ Tip


The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on and continue tracing along the gumline.

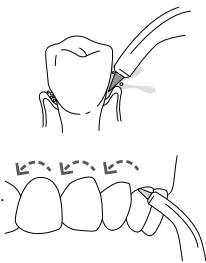


### Tongue Cleaner Tip

To use the tongue cleaner tip, set to lowest pressure setting, place tip in the center/middle of your tongue about half way back and start the Water Flosser. Pull forward with light pressure. Increase pressure as you prefer.

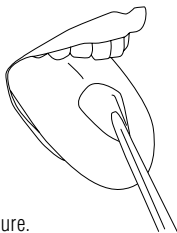

### Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

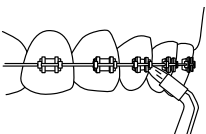

### Toothbrush Tip

Place the toothbrush tip in mouth with the brush head on the tooth at the gum line. The toothbrush tip can be used with or without toothpaste. Turn the Water Flosser on so water flows through the tip. Using a **light pressure** (bristles should not bend) massage the brush back and forth with very short strokes – much like you would with a manual toothbrush.



### Plaque Seeker® Tip

To use the Plaque Seeker® tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.



### Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

## Cleaning

Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

The water reservoir is removable for easy cleaning and is top rack dishwasher safe.

**NOTE:** Remove reservoir valve before placing reservoir in dishwasher.



Remove



Replace

## Removing Hard Water Deposits

Add 2 tablespoons of white vinegar to a full reservoir of warm water. Set pressure control dial to high, point the handle and tip into sink. Turn unit ON (I) and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

## Leakage

If leaking due to a damaged hose, contact Waterpik International, Inc. UK Customer Service Centre at +44 (0) 844 335 6354 for a self service repair kit.

If leakage occurs from the reservoir, ensure the valve is right side up or contact Waterpik International, Inc. UK Customer Service Centre for a replacement valve.



## Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your service needs please call the Waterpik International, Inc. UK Customer Service Centre at +44 (0) 844 335 6354. Some replacement parts and accessories for the Water Flosser are also available at [www.waterpik.co.uk](http://www.waterpik.co.uk).

For questions or problems please call the customer service centre before returning the product. Please do not send your unit to the Waterpik International, Inc. address as this can cause lengthy delays.

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Still have questions? Go to [www.waterpik.co.uk](http://www.waterpik.co.uk).



Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for three years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.







**Water Pik, Inc.**

1730 East Prospect Road  
Fort Collins, CO 80553-0001  
USA

**Waterpik, International, Inc.**

Castle Court, 41 London Road  
Reigate, Surrey RH2 9RJ  
United Kingdom

**[www.waterpik.co.uk](http://www.waterpik.co.uk)**

Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Benelux, Canada, Chile, China, Hong Kong, India, Japan, Korea, Mexico, New Zealand, Norway, Pakistan, Russian Federation, Singapore, South Africa, Switzerland, Taiwan, Ukraine, and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Australia, Brazil, Canada, China, EU, Hong Kong, Japan, Mexico, Korea, Russian Federation, Switzerland, Trinidad & Tobago, and the United States.

Pik Pocket™ is a trademark of Water Pik, Inc.

Plaque Seeker® is a trademark of Water Pik, Inc. registered in Canada, China, EU, and the United States.

Designed in USA. Made in China.

